

PRODUCT INFORMATION

IAS3119

2019EC

Please note that values for Plasma Renin Activity are not available for
Immunoassay Speciality II Level 3 lot 2019EC.

CCS6875

IMMUNOASSAY SPECIALITY II - LEVEL 3 (IA SPECIALITY II LEV 3)

CAT. NO. IAS 3119

LOT NO. 2019EC

SIZE: 5 x 1 ml

EXPIRY: 2022-10-28

GTIN: 05055273207347

INTENDED USE

This product is intended for *in vitro* diagnostic use, in the quality control of the accuracy of Immunoassays on clinical chemistry systems. This material can be used to monitor the control of accuracy or the control of reproducibility of immunoassays.

SAFETY PRECAUTIONS AND WARNINGS

For *in vitro* diagnostic use only. Do not pipette by mouth. Exercise the normal precautions required for handling laboratory reagents.

Human source material, from which this product has been derived, has been tested at donor level for the Human Immunodeficiency Virus (HIV 1, HIV 2) antibody, Hepatitis B Surface Antigen (HbsAg), and Hepatitis C Virus (HCV) antibody and found to be NON-REACTIVE. FDA approved methods have been used to conduct these tests.

However, since no method can offer complete assurance as to the absence of infectious agents, this material and all patient samples should be handled as though capable of transmitting infectious diseases and disposed of accordingly.

Health and Safety Data Sheets are available on request.

STORAGE AND STABILITY

OPENED: Store refrigerated (+2°C to +8°C). In reconstituted serum Renin is stable for 5 days, Procalcitonin is stable for 1 day, Gastrin and Calcitonin are stable for 8 hours at +2°C to +8°C if kept capped in original container and free from contamination. The control is stable if frozen once for 28 days at -20°C, except for the Siemens Advia Centaur CALCT assay (Calcitonin) which is stable for 14 days at -20°C. Only the required amount of product should be removed. After use, any residual product should NOT BE RETURNED to the original vial.

UNOPENED: Store refrigerated (+2°C to +8°C). Stable to expiration date printed on individual vials.

Bacterial contamination of the reconstituted serum will cause reductions in the stability of many components. If bacterial contamination is suspected, the vial should be discarded and a fresh vial reconstituted.

PREPARATION

Immunoassay Speciality II is supplied lyophilised.

1. Carefully reconstitute each vial of lyophilised serum with exactly 1 ml of distilled water at +15°C to +25°C. Close the bottle and allow to stand for 30 minutes before use. Ensure contents are completely dissolved by swirling gently. Avoid formation of foam. Do not shake.
2. Refer to the Control section of the individual analyser application.
3. Refrigerate any unused material. Prior to reuse, mix contents thoroughly.

MATERIALS PROVIDED

Immunoassay Speciality II - Level 3 5 x 1 ml


MATERIAL REQUIRED BUT NOT PROVIDED

Volumetric pipette

VALUE ASSIGNMENT

Each batch of Immunoassay Speciality II is submitted to a number of reference laboratories and values are assigned from a consensus of results obtained by these laboratories, using a unique statistical analysis. With each batch, a control range is provided for individual parameters and each parameter method. The control range is equivalent to the assigned mean ± 2 S.D. This results in extremely accurate values, which may be confidently used by laboratories to ensure the accuracy of their methods.

If a method is unavailable, contact Randox Laboratories - Technical Services, Northern Ireland, tel: +44 (0) 28 9445 1070 or email Technical.Services@randox.com.

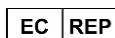

Randox Teoranta, Meenmore,
Dungloe, Donegal,
F94 TV06, Ireland

03 Mar '21 me

IMMUNOASSAY SPECIALITY II - LEVEL 3 (IA SPECIALITY II LEV 3)

Cat. No. IAS3119 Lot. No. 2019EC Size 5 x 1ml Expiry 2022-10-28

| Range | | | | | |
|---------------|--------|--------|------|------|-------------------------------|
| Analyte | unit | Target | low | high | methods |
| Calcitonin | pmol/l | 96.5 | 72.4 | 121 | Siemens Immulite 2000/2500 |
| | pg/ml | 329 | 247 | 411 | |
| | pmol/l | 148 | 111 | 185 | Roche Elecsys/Cobas/Modular |
| | pg/ml | 505 | 379 | 631 | |
| Gastrin | pmol/l | 175 | 131 | 219 | Siemens Immulite 2000/2500 |
| | pg/ml | 365 | 273 | 457 | |
| Procalcitonin | µg/l | 20.2 | 15.2 | 25.3 | BioMerieux Vidas |
| | µg/l | 16.6 | 12.5 | 20.8 | Roche Elecsys/Cobas/Modular |
| | µg/l | 16.8 | 12.6 | 21.0 | Abbott Architect Brahms PCT |
| | µg/l | 7.92 | 5.94 | 9.90 | Shanghai I-Reader |
| Renin | mIU/l | 365 | 274 | 456 | Diasorin Liaison direct Renin |
| | pg/ml | 215 | 161 | 269 | |